

Xiaomi Electric Toothbrush T700

User Manual



Safety Instructions

When using the toothbrush for the first time, it may cause itch or minor gum bleeds. This is normal when the gums are not accustomed to the efficient cleaning of an electric toothbrush. It is recommended to use custom mode when your gums are not used to it yet. Normally, gum bleeds will stop within one to two weeks. Stop using the toothbrush and consult a dentist if it persists.

WARNING

- Users undergoing orthodontics, dental repair, or dental implant should use the product under the guidance of their dentists.
- This toothbrush complies with the safety standards for electromagnetic devices. If you have any questions about using the device alongside pacemakers or other implanted devices, please consult a doctor or the manufacturer of the implanted device before use.
- This toothbrush is designed to clean teeth, gums, and tongue only. Do not use it for any other purposes. If you feel any discomfort or pain, stop using the toothbrush and consult your dentist.
- Do not immerse the charging dock in water or any other liquid.
- Do not plug the charging dock into an outlet when there is liquid on the plug.
- This toothbrush contains a built-in rechargeable battery, and should never be exposed to fire, nor charged, used, or stored in a high-temperature environment.
- Do not hit or drop the toothbrush. Stop using it if it has been damaged due to impact or dropping.
- Do not attempt to modify, disassemble, or repair the toothbrush by yourself.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- This toothbrush is equipped with a built-in battery that cannot be removed or replaced. Do not disassemble or modify the battery by yourself.

- A minimum separation distance of 0 cm must be maintained between the user's body and the device to comply with the RF exposure requirements in Europe.
- Only use the charger provided with the product. If supplied, only use the original USB power supply unit.
- Do not use the charger outdoors or near heated surfaces.
- For the purposes of recharging the battery, if supplied, only use the detachable supply unit provided with this appliance: A520-050100U-EU1.
- Keep the product and the battery away from fire and do not expose them to direct sunlight. Do not expose the appliance or battery to excessive temperatures.
- Do not open, modify, pierce, damage or dismantle the product or battery to prevent battery from heating up or releasing toxic or hazardous substances. Do not short-circuit, overcharge or reverse charge the battery.
- If the battery is damaged or leaking, avoid contact with the skin or eyes. If this occurs, immediately rinse well with water and seek medical care.
- Do not place product and its battery in microwave ovens or on induction cookers.
- This appliance contains a battery that is non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.

CAUTION

- Check the toothbrush for obvious damage before using it. Do not use it if it is severely damaged.

- Before plugging into a power outlet, make sure the voltage of the local power outlet matches the rated voltage of the product.

- Do not use the charging dock outdoors in direct sunlight.
- Do not use the toothbrush while bathing or showering.
- Do not share a toothbrush head with others.
- Do not wrap the charging cable around the charging dock while it is plugged into an outlet.


- Do not place any metal objects on the charging dock while it is plugged into an outlet.

Note: Illustrations of product, accessories, and user interface in the user manual are for reference purposes only. Actual product and functions may vary due to product enhancements.

For detailed e-manual,
please go to www.mi.com/global/service/userguide

Product Overview

Read this manual carefully before use, and retain it for future reference.

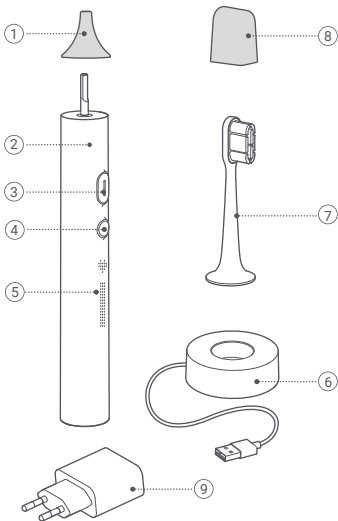


- ① Shaft Protection Cover
- ② Toothbrush Handle
- ③ On/Off Button
- ④ Mode/Power Button

- ⑤ Smart LED Screen
- ⑥ Charging Dock
- ⑦ Toothbrush Head
- ⑧ Toothbrush Head Cover
- ⑨ Power Adapter

Installing the Toothbrush Head

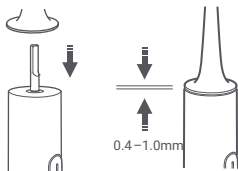


Make sure that the toothbrush head is installed in the same direction as the front of the toothbrush handle, and make sure there is a slight gap in between the toothbrush head and the handle. (The gap is designed to prevent the toothbrush head from colliding with the handle during high-frequency vibration, so as to generate more power from the motor.)



How to Use

	Turning on/off
	On/Off button
	Turning on: Press the button to turn on the toothbrush. The toothbrush will enter Standard Mode by default when powered on for the first time.
	Turning off: Press the button to turn off the toothbrush.
	Changing mode/power
	Mode/Power button
	Changing modes: When the toothbrush is powered on, press the button to change modes.
	Stepless power adjustment: When the toothbrush is powered on, press and hold the button to adjust the power in a stepless manner.
	Notes: The brushing timer restarts every time the power is changed. Every power adjustment is only made to the current mode. Power in other modes is not affected.

Introduction to brushing modes

-  Standard Mode: Suitable for daily oral hygiene.
-  Gentle Mode: Suitable for sensitive teeth.
-  Custom Mode: A two-minute advanced mode by default. More features and modes can be configured in the Mi Home/Xiaomi Home app including brushing duration, splash prevention for ten seconds, beginner mode, sensitive mode, advanced mode, additional features, etc. These settings are available in Custom Mode.



Restoring factory settings

Press and hold the on/off button and the mode/power button at the same time for five seconds, and the smart LED screen will turn white until the buttons are released. This will restore the product to the factory settings.

Note: To fully experience functions of the toothbrush's app, please connect to the Mi Home/Xiaomi Home app in time after restoring factory settings.

Healthy brushing

- Position the toothbrush head at an angle of 45° to the gum line, and slowly move it along the teeth and gum line.
- While brushing, the toothbrush pauses briefly every 30 seconds to prompt the user to change the brushing position.
- When there is no active pause during use, the toothbrush will automatically turn off once the brushing mode ends. The brushing duration can be set via the Mi Home/Xiaomi Home app.

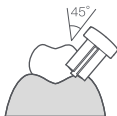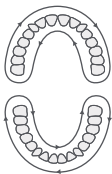

Screen icons



Standard mode



Gentle mode



Custom mode



Power indicator



Preparing Bluetooth connection
(Indicators will breathe.)



Connecting via Bluetooth
(Indicators will sequentially turn on and off.)



Bluetooth connected successfully



Brushing score



Brushing score is over 95



Brushed for less than one minute

- Power indicator: The stronger the power, the more the indicator lights up.
- Brushing score: After every brushing session is completed, the screen will display the score for the session. If the brushing score is over 95, then "95" will be displayed. If the brushing session lasts for less than one minute, then "..." will be displayed. The Mi Home/Xiaomi Home app does not record the brushing data for sessions less than one minute. It is recommended to brush for at least one minute every time.

Screen-on notifications

- When the toothbrush is powered on for the first time, the screen-on notification feature will be enabled automatically. Afterward, whenever you pick up the toothbrush, the screen will automatically light up and display the current brushing mode.
- If a nickname or anniversary date has been set in the Mi Home/Xiaomi Home app, then the screen will display the nickname or the anniversary notification on the set date when the toothbrush is picked up.

Notes:

- The brushing mode will not be displayed when the nickname or anniversary notification is displayed.
- When the toothbrush is carried during travel, the screen-on notification feature automatically will be off temporarily if continued bumping is detected. It will be on again when the toothbrush is turned on again or there is no movement for 30 minutes.


Cleaning the toothbrush body

- Clean the toothbrush body with a dampened soft cloth at least once a week.
- Toothpaste should never get on the body of the toothbrush; if it does, immediately wash the entire body of the toothbrush under a gentle stream of running water at room temperature and then dry the entire body of the toothbrush with a soft cloth.

Cleaning the toothbrush head

- After brushing, rinse the toothbrush head with running water.
- After rinsing, put on the toothbrush head with the protective cover to keep it clean and hygienic.


Note: It is recommended to replace the toothbrush head every three months.

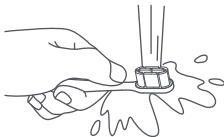

Smart Functions

Overpressure notification

The toothbrush is equipped with an overpressure protection system that detects pressure in real-time. When too much pressure is applied during teeth brushing, the overpressure notification indicator will blink in red quickly and the toothbrush will vibrate alternately in different amplitudes to notify the user.

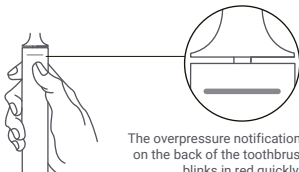

The overpressure notification indicator on the back of the toothbrush handle blinks in red quickly.

Connecting with Mi Home/Xiaomi Home app










- This product works with Mi Home/Xiaomi Home app*. Control your device with Mi Home/Xiaomi Home app.
- Scan the QR code to download and install the app. You will be directed to the connection setup page if the app is installed already. Or search "Mi Home/Xiaomi Home" in the app store to download and install it.
- Open Mi Home/Xiaomi Home app, tap "+" on the upper right, and then follow prompts to add your device.



BF149176


Note: The version of the app might have been updated, please follow the instructions based on the current app version.

Regular Charging

Toothbrush Status	Battery Level Indicator Status	Battery Level
Using	 Quickly blinking	<25%
	 Quickly blinking (unable to turn on)	Battery depleted
Charging	 Fully lighted (charging completed)	100%
	 The 5th bar slowly blinking	100%–80%
	 The 4th bar slowly blinking	80%–60%
	 The 3rd bar slowly blinking	60%–40%
	 The 2nd bar slowly blinking	40%–20%
	 The 1st bar slowly blinking	<20%
	 Slowly blinking	0%

- The toothbrush charging dock is designed for the T700 and is not compatible with other products.
- When the toothbrush cannot be turned on after the battery is depleted, it will take approximately six hours to fully charge it. The time required for fully charging will depend on the battery level when charging starts.

*If the toothbrush battery has been depleted and was not charged for a long period of time, the charging icon may blink slowly the next time it is charged. This is normal.

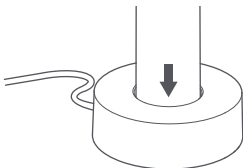

Specifications

Name: Xiaomi Electric Toothbrush T700

Model: MES604

Rated Voltage: 3.7 V ---

Rated Power: 4.2 W

Rated Input : 5 V --- 1 A

Charging Time: Approx. 6 hours

IP Rating: IPX7

Wireless Connectivity: Bluetooth 5.0 2.4 GHz

Operation Frequency: 2402–2480 MHz

Maximum Output Power: <4 dBm

Power Adapter

Information Published	Value and Precision	Unit
Manufacturer	Jiangxi Jian Aohai Technology Co., Ltd.	-
Model Identifier	A520-050100U-EU1	-
Input Voltage	100–240	V \sim
Input AC Frequency	60/50	Hz
Output Voltage	5.0	V ---
Output Current	1.0	A
Output Power	5.0	W
Average Active Efficiency	73.62	%
Efficiency at Low Load (10%)	64.59	%
No-load Power Consumption	0.075	W

Troubleshooting


Issues	Possible Causes	Solutions
Unable to turn on	Charging	Remove the toothbrush from the charging dock
	Battery level is too low	Plug into the charging dock to charge the toothbrush
Unable to install the toothbrush head onto the shaft	There is something stuck inside the toothbrush head	Use a tool to remove the object from the toothbrush head
	There is something stuck on the shaft	Clean the shaft
	Not a genuine toothbrush head	Replace it with a genuine toothbrush head
Excessive power consumption and charging frequently	Connected to the phone via Bluetooth for a long time	Turn off Bluetooth when the toothbrush is not used
	Used too frequently	Plug into the charging dock to charge the toothbrush
Toothbrush head vibration is weak	The toothbrush head is not installed properly	Reinstall the toothbrush head
	There is something stuck between the toothbrush head and the handle	Remove the toothbrush head and remove the object
Cannot connect to the Mi Home/Xiaomi Home app	Bluetooth is not enabled on your phone	Enable Bluetooth on your phone and control it via the Mi Home/Xiaomi Home app
	Toothbrush is too far away from the phone	Move the toothbrush close to the phone
	Toothbrush is already turned on	Turn off the toothbrush and try connecting again
	Battery level is too low	Try connecting again after charging

Environmental Protection

When the toothbrush reaches the end of its service life, please remove, dispose of, or recycle the built-in battery in accordance with local laws and regulations of the country or region it is used in. If the battery cannot be removed, please take the toothbrush to an official designated recycling center.

CAUTION

The toothbrush must be turned off before removing the battery. Before removing the battery, make sure it is depleted. You can make the toothbrush remain on until it stops working, and repeat this process until the battery runs out of power.

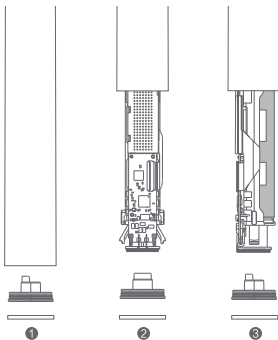

How to remove the battery

- ① Remove the decoration and the cover at the bottom of the toothbrush handle.
- ② Remove the inner holder.
- ③ Cut off the cable connecting the battery and remove the battery.

Note: Do not replace the battery on your own. Removing the battery forcibly may damage the function of the toothbrush, please be cautious.

Regulatory Compliance Information

WEEE Disposal and Recycling Information



All products bearing this symbol are waste electrical and electronic equipment (WEEE as in directive 2012/19/EU) which should not be mixed with unsorted household waste.

Instead, you should protect human health and the environment by handing over your waste equipment to a designated collection point for the recycling of waste electrical and electronic equipment, appointed by the government or local authorities. Correct disposal and recycling will help prevent potential negative consequences to the environment and human health. Please contact the installer or local authorities for more information about the location as well as terms and conditions of such collection points.

EU Declaration of Conformity



Hereby, Wuxi Qinghe Xiaobei Technology Co., Ltd. declares that the radio equipment type MES604 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address:
<http://www.mi.com/global/service/support/declaration.html>

Manufactured for: Xiaomi Communications Co., Ltd.

Manufactured by: Wuxi Qinghe Xiaobei Technology Co., Ltd.
(a Mi Ecosystem company)

Address: Room 1106, 11th Floor, NO.2 of Building 530, Binhu District,
Wuxi, China

For further information, please go to www.mi.com

User Manual Version: V1.0

Importer:

Beryko s.r.o.

Na Roudné 1162/76, 301 00 Plzeň

www.beryko.cz